



RC-SL series

Minor surgery light

1. Double domes both are 160000Lux with 5 inch LCD touch control panel
2. Has 10 steps adjustable illumination
3. Has 5 steps adjustable color temperature
4. Has Endo mode
5. Deep illumination
6. PWM dimming technology, no flickering
7. LED cold light, head temperature rise less than 1°C
8. High color rendering index ensures clear views of the tissue and Organ
9. Removable sterilizer handle
10. Light field adjusting
11. Beautiful ultra-thin lamp head design
12. Has R9 Red light mode adjustment, red light mode is a color rendering index specifically for red, can accurately see the subtle differences in red tissue; And Green light mode is suitable for cardiothoracic surgery, neurosurgery and other surgical operations requiring high color rendering index
13. Illuminance can be direct adjusted in 25% and 60%



OPTIONAL MODEL

**Model No. RC-SL700M Touch
RC-SL500M Touch**

Mobile type
(Backup battery is optional)

RC-SL700+700 touch






RC-SL700+700 touch



RC-SL700+500 touch

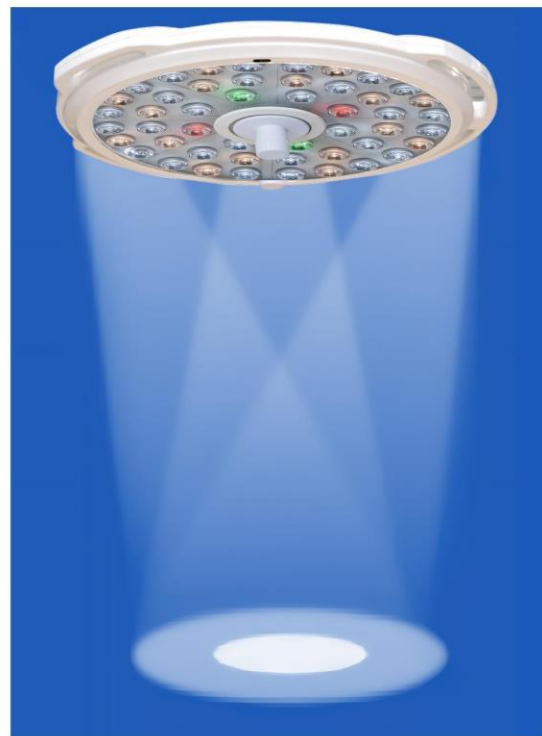

Has R9 Red light mode adjustment, red light mode is a color rendering index specifically for red, can accurately see the subtle differences in red tissue; And Green light mode is suitable for cardiothoracic surgery, neurosurgery and other surgical operations requiring high color rendering index

Has 10 steps adjustable illumination
Has 5 steps adjustable color temperature
Has Endo mode
Deep illumination



RC-SL700+500+TV touch

Double domes main dome is 160000Lux and satellite dome is 140000Lux with 5 inch LCD color touch control panel
Build-in 3 million pixel Sony HD Camera, 1080P resolution 24inch Dell monitor, large capacity video recorder
Optional electric focus adjustment function
Optional illuminance compensation function

Main Technical Data

RC-SL700+700touch	Minor Surgery LED Light
Model	RC-SL700+700touch
Head Diameter	700mm/700mm
No. Of LED Blubs	80pcs+80pcs
Illumination	160000Lux/160000Lux
CRI	Ra>95
Color Temperature	3500-5000K
Illumination Depth	1200mm
Spot Diameter	120-300mm
LED Life	>50,000 Hours
LED Brand	OSRAM GWCS8PM1.CM
UPS	Optional



RC-SL700 Touch

Ceiling type single dome 160000Lux with 5 inch LCD touch control panel